BEOSIGNER®

PRODUCT MANUAL

INDEX

Deutsch	2
English	9
Français	16
Magyar	23
Nederlands	30
Italiano	37
Español	44
Polski	51
اللغة العربية	59

VITARIGHTS
WISSEN SCHAFFT DIALOG

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BEOSIGNER® INSTRUCTION MANUAL

Congratulations on your new Beosigner®

Please read this instruction manual and the safety instructions for the Beosigner® carefully before use and keep them for future reference.

Warning! Direct eye contact should be avoided. If the laser beam hits the eye, close your eyes immediately and turn your head away from the beam. Do not use optical instruments such as a magnifying glass to view the radiation source.

The Beosigner® is not a toy. Parents should keep the Beosigner® out of reach of children. If parents allow their children to use the Beosigner®, they should carefully supervise their handling of the device.


This product can be used by persons with reduced physical, sensory or mental capabilities if they have been given supervision or instruction concerning the use of the product in a safe way and understand the hazards involved.

Attention! Please keep the Beosigner® away from microwaves, as microwaves can impair its functioning. Do not expose the Beosigner® to steam or moisture.

Never open the Beosigner® without authorization. Modifications and repairs to the device may only be carried out by the manufacturer.

Treatment with the Beosigner® is very cell-compatible, as scientific studies* have shown. However, as there is insufficient experience to date in connection with use on the body in pregnant women and infants (up to 2 years of age) and in the case of cancer and organ transplants, it cannot be recommended at present. Water, food and drinks can always be vitalized with the Beosigner®.

**Cf. Wiegand C, Dirksen A, Tittelbach J. Treatment with a red-laser-based wound therapy device exerts positive effects in models of delayed keratinocyte and fibroblast wound healing. Photodermatol Photoimmunol Photomed. 2023;00:1-14. doi:10.1111/phpp.12926.*

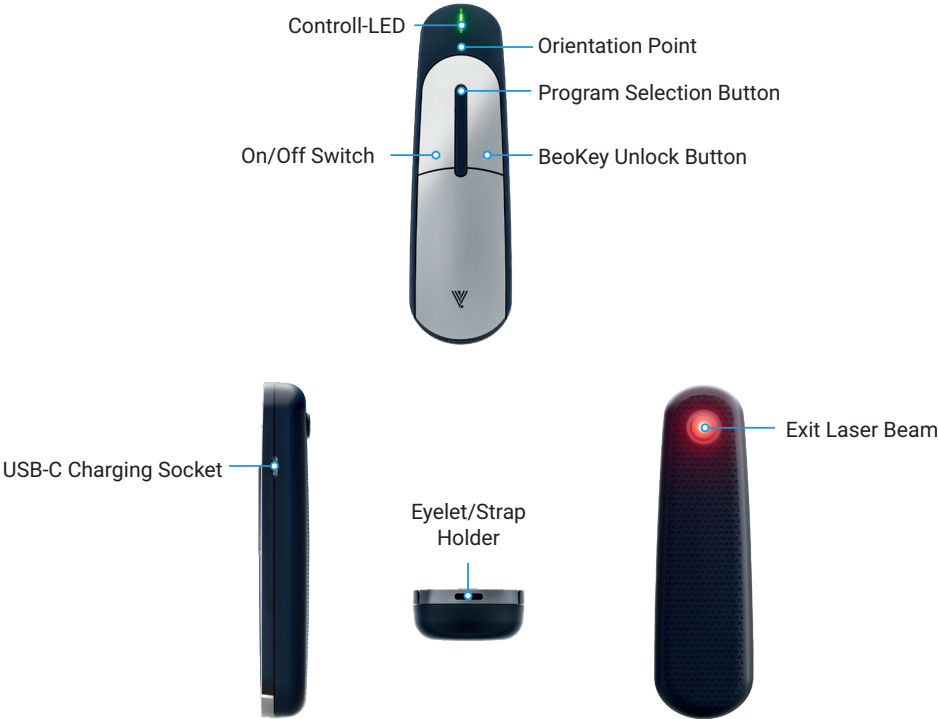Charging

Please fully charge your Beosigner® before first use. (Charging time approx. 10 hours).

To charge, connect the USB-C cable of the charger to the charging socket of your Beosigner® and the charger to a power outlet. The control LED flashes red during charging and goes out as soon as the Beosigner® is fully charged. The Beosigner® cannot be switched on during the charging process. Do not store the Beosigner® for a longer period of time (maximum 3 months) without charging the device.

Switching on and off

Switching on: Press the on/off switch (left button of the Beosigner®).

The Beosigner® is now in standby mode, the control LED lights up green.

Immediately after switching on, you can recognize the battery charge level by the color of the Control LED:

Control LED at power on	Battery charge level
permanently green	full or still sufficient (over 20 %)
flashes green and blue alternately for 2 sec	below 20 %
flashes green and red alternately for 2 sec	below 10 %

Switching off: Regardless of which program you are in, the Beosigner® can be switched off at any time. To do this, press the on/off switch again (left button of the Beosigner®).

Note: After 1 minute in standby mode, the Beosigner® switches off automatically.


The Beosigner® has 2 different programs

Program 1: Normal function (10 minutes)

Switch on the Beosigner® (press on/off switch).

Briefly press the program selection key (middle key). The Control LED lights up red and the laser light on the underside of the Beosigner® is switched on.


The program runs for a total of 10 minutes. After the following time intervals, the control LED flashes briefly and a short acoustic signal is heard: 15 sec – 30 sec – 45 sec – 60 sec – then every 60 sec.

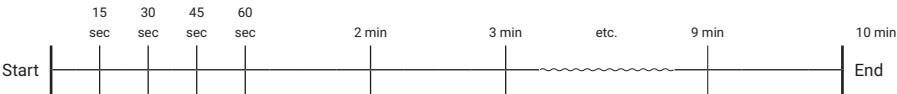Program 2: Interval switching (total 10 hours)

Switch on the Beosigner® (press on/off switch).

Press the program selection key (middle key) for approx. 3 seconds until the control LED flashes blue. Release the program selection key. The control LED lights up blue and the laser light on the underside of the Beosigner® is switched on for 10 minutes. Then the Beosigner® switches off for 50 minutes (sleep mode, control LED pulses blue), then on again for 10 minutes. This is repeated for the entire runtime of 10 hours.

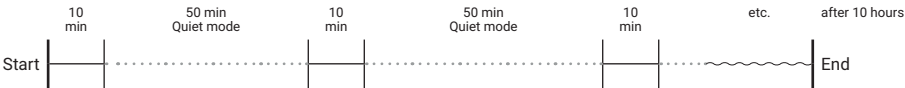Program 1 can be ended prematurely at any time by briefly pressing the program selection key (middle key). After the program has finished, the laser light goes out and the Beosigner® switches to standby mode.

Program 2 can be ended at any time by pressing the on/off key.

Switching the advisory tones on and off


Switch on the Beosigner® (press on/off switch). The Beosigner® is now in standby mode.

To permanently switch the advisory tones on or off, press and hold the program selection key (middle key) until the control LED first flashes blue and then green after a few more seconds. Release the program selection key.

Orientation point

The Beosigner® features an orientation point: a small elevation that can be easily felt with the finger, which makes it easier to position the laser light that emerges exactly opposite to it on parts of the body that cannot be seen.

Orientation Point

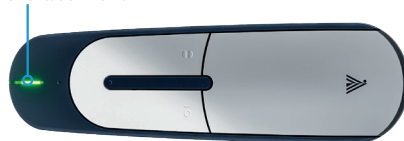Eyelet/strap holder

The supplied ribbon can be threaded into the eyelet and used in conjunction with a ribbon or chain of your choice to hang the Beosigner® around your neck.


Eyelet/Strap Holder

Cleaning and care


Do not clean the Beosigner® while it is switched on or connected to the charging cable. Clean the Beosigner® with a soft, dry cloth. Remove stubborn dirt with a slightly damp cloth. A soft cotton swab or similar can also be used to clean the lens (laser beam exit); however, avoid exerting strong pressure on the lens.

Remove the protective film from the BeoKey

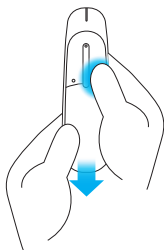Releasing the BeoKeys

Press and hold the unlock button (right button of the Beosigner®). The BeoKey can now be easily slid out of the Beosigner® horizontally without using force.

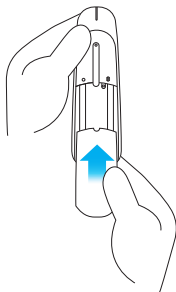Inserting the BeoKeys

Insert the BeoKey at the end of the guide rail and carefully push it in horizontally – without applying force – until it engages in the Beosigner®. Avoid tilting at all costs, both laterally and vertically. The force applied could irreparably damage the BeoKey and the Beosigner®.

Safety instructions for the BeoKeys

We recommend that you do not wear the BeoKeys directly on the skin for long periods of time, as contact with creams or other cosmetics, for example, could cause optical changes to the surface.

Manipulation of the retaining bolt, opening the Beosigner® housing or forcible removal of the BeoKey using tools will invalidate any warranty claims.

Disposal

The symbol “crossed-out dustbin” indicates that the device or the battery must not be disposed of with household waste, but must be taken to specially set up collection points, recycling centers or disposal companies.

Trademark notice

BEOSIGNER® is a registered trademark of Vitarights Innovations GmbH.

The product complies with the standard DIN EN 60825-1:2022-07.

Please keep this instruction manual with the packaging and make them available to anyone using the Beosigner®.

In Switzerland, the use of class 2 laser devices is prohibited by law. This also applies to the import and transit, offering, sale and use – and therefore also to transit through Switzerland (e.g. by plane, train, motor vehicle, bicycle, on foot) with a Beosigner® class 2 laser.

Technical Data

- Laser class 2
- Maximum power of laser radiation with $P \leq 1 \text{ mW}$
- Wavelength: 640 nm - 660 nm

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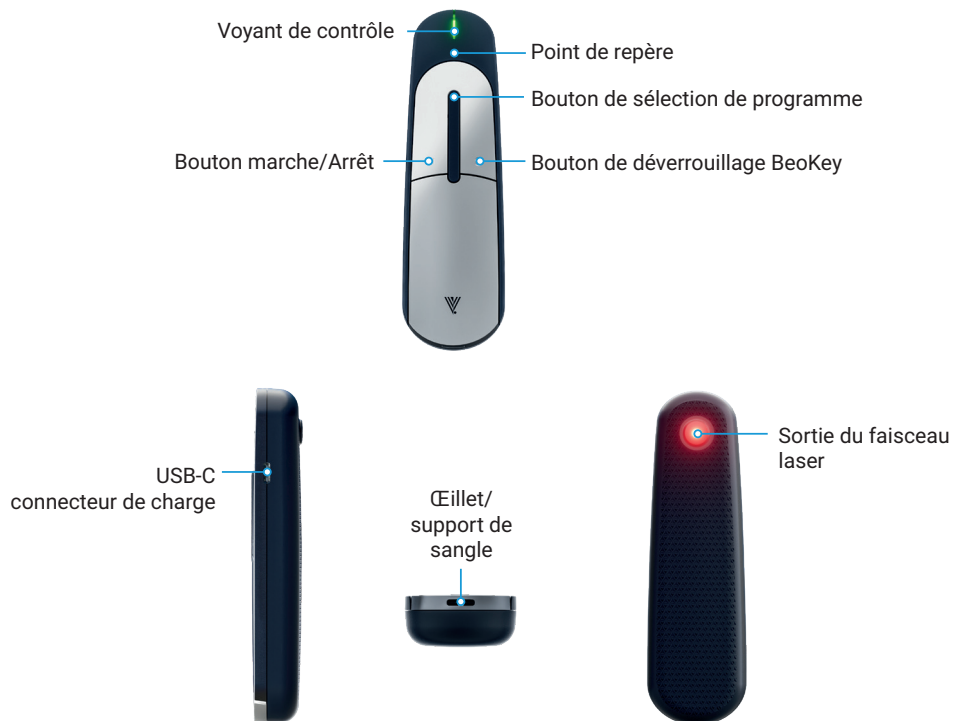
*Cf. Wiegand C, Dirksen A, Tittelbach J. Treatment with a red-laser-based wound therapy device exerts positive effects in models of delayed keratinocyte and fibro-blast wound healing. Photodermatol Photoimmunol Photomed. 2023;00:1-14. doi:10.1111/phpp.12926.

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